RED

AESTHETIC

POLYMER/MONOMER

INSTRUCTIONS FOR USE

GEBRAUCHSINFORMATION

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HEAT-CURING DENTURE BASE MATERIAL

Dear customer,

The product that you have purchased qualifies as a medical device under European Directive 93/42/EEC. Please ensure that you file the lot number and product name in your administrative system on receipt of the product. For each piece of work that you produce, please make a note of all the materials used along with all the appropriate lot numbers. For your own protection as well as for the protection of your patients and the environment, please also follow the important guidelines in the following instructions for use.

Product description

AESTHETIC RED sets new quality standards in acrylic heat-curing denture base materials where processing, color and dimensional stability are concerned. The dental technician has different opaque shades, transparent pink and clear to choose from. Various polymerisation methods are open to the user, permitting a large number of application possibilities.

Composition / Powder

Polymethyl methacrylate, softening agent, benzol peroxide, pigments

Composition / Liquid

Methyl methacrylate, dimethacrylate, catalyst

Storage instructions

Store the material in a cool, dark, well-ventilated place.

Storage temperature: 2–28 °C / 36 –82 °F. Do not use the material after the expiry date.

Keep out of the reach of children.

Indication

- Full dentures
- Partial dentures
- Combination dentures
- Relining

Contraindication

Avoid direct contact with unpolymerised material within the oral cavity. If a patient is known to be allergic to any of the ingredients in AESTHETIC RED, the material must not be used.

Side effects

No systemic side effects have been reported to date. In individual cases, local allergic reactions to PMMA/MMA-based denture base materials have been reported.

Danger warnings

- Monomer contains methyl methacrylate (MMA)
- MMA is an irritant and easily flammable (flash point: + 10 °C / 50 °F)
- MMA and its vapours are irritating to the eyes, skin and respiratory system
- May cause sensitisation upon skin contact
- Do not inhale vapours
- Keep away from sources of ignition – no smoking
- Do not empty into drains

- Avoid contact of the skin with monomer and uncured material. Many commercial gloves, e.g. those made of latex or vinyl, are not monomer-resistant and therefore do not provide protection against the sensitising effect of methacrylates.
- Wear a mask when grinding and use a suction removal system.
- Always use a spatula when handling the mixture.
- Take precautionary measures against static discharges.

Warning

This material has been developed solely for use in dentistry and must be handled strictly in accordance with the instructions for use. The manufacturer cannot be held liable for any damages arising as a result of failure to observe the instructions for use or the stated area of application. The same applies in the event that the product is mixed or processed with other manufacturers' products. The user shall be solely responsible for testing the material with respect to its suitability prior to use for any purpose other than those explicitly stated in the instructions.

Safety data sheets can be found on our website at www.candulor.com.


INSTRUCTIONS FOR USE

1. Preparation

The boiled out precasts must be warm (approx. 40–50 °C / 104 –122 °F) and the plaster surfaces thoroughly water-saturated. Isolate the surfaces with two applications of ISO-K and allow to dry thoroughly. Roughen the teeth and wet with monomer to ensure proper bonding with the denture base.

2. Dosage

Ideal mixing ratio for one denture:
22 g polymer : 10 ml monomer

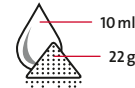The dosage system ensures an ideal mixing ratio and minimum polymerisation shrinkage of AESTHETIC RED. The scale on the monomer cylinder is in millilitres, on the polymer beaker in grams.

Too high a monomer content will alter the shade, consistency and handling properties of the material. This may also inhibit polymerisation, possibly leading to irritation of the mucous membranes and to porosity.

3. Mixing

Mix together the powder and liquid based on the recommended ratio or as required, using a spatula, and blend thoroughly for approx. 20 sec. Cover the mixing beaker and leave to prove for approx. 8–10 minutes at room temperature (23 °C / 73 °F).

The powder and liquid must be mixed together thoroughly for approx. 20 seconds using a hard metal spatula.

4. Working time

After the proving time, as soon as the material is no longer sticky, it can be processed for approx. 20 minutes at 23 °C / 73 °F.

The proving and working time is variable, depending on the mixing ratio selected!

5. Pressing

Place an ample amount of the mixture in one half of the flask which must be at body temperature and which you have previously wetted and isolated with ISO-K.

Close the flask carefully, place it under the press and apply a pressure of 80 bar.

Do not change the pressure intensity!

6. Polymerisation

The dental technician can choose between different methods of polymerisation:

Standard processing

Place the closed flask in cold water, heat to 100 °C (212 °F) and boil for 45 minutes.

Other methods of processing

- Place flask in cold water, heat to 70 °C (158 °F) and leave it at this temperature for 30 minutes. Then heat to 100 °C (212 °F) and boil for 30 minutes.
- Place flask in hot water (70 °C / 158 °F) and leave it at this temperature for 60 minutes. Then heat to 100 °C (212 °F) and boil for 30 minutes.
- Place flask directly into boiling water, bring the water back to the boil and boil for 40 minutes. Only suitable for small to medium-sized dentures!
- Place flask in cold water, heat to 80 °C (176 °F) and polymerise for 10 hours. Then turn off the heat source and leave flask to cool in the same water overnight.
- Leave flask to polymerise in a drying cabinet at 80 °C (176 °F) for 10 hours.

By using the ideal mixing ratio and the recommended methods, a residual monomer content of < 2.2% can be achieved.

7. Cooling

Leave the flask to stand in air for at least 30 minutes, then cool completely in cold water. Allow the flask to cool completely before opening.

8. Removal and finishing

When the flask has cooled completely, open it and remove the plaster. Check the occlusion of the denture, then finish and polish in the usual way.

Repairs and corrections of AESTHETIC RED

Repairs and corrections may be carried out using AESTHETIC BLUE (cold-curing denture base material). The contact surfaces must be well roughened and wetted with monomer.

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🏢 Candulor AG